ANTICELL MASSAGER

Instructions manual • Instrukcja obsługi • Návod k použití • Installieremise juhend • Įdiegimo instrukcija
• Installasjonsmanual • Руководство по установке • Asennusohje • Uzstādīšanas instrukcija
• Εγχειρίδιο χρήσης • Gebrauchsanleitung • Használati útmutató • Manual de instruções

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Cellulite is an abnormal distribution of fat tissue leading to changes in the subcutaneous tissue. The degree of skin changes may be visible to the naked eye or only noticeable after squeezing a skin fold. Uneven skin texture can be effectively removed with the right massage.

Device

- ① Massage head
- ② LED display of operating modes and battery charge level
- ③ ON/OFF touch button – change of operating modes
- ④ USB-C charging port

Massage pads

1. Mesh cover - can be used as an addition to massage attachments, preventively, as protection against hair or other small materials getting entangled in the massager mechanism.
2. Flat head - deep massage, head ideal for relaxing massage of the whole body
3. Head with lipoactive balls - the action of this attachment promotes the movement of the head on the surface of the body and generates a feeling of continuity and fluidity, creating an effective lipoactive massage.
4. Finger head - simulation of finger massage, for the whole body
5. Wave head - for massage of soft tissues, relaxes sensitive muscles, creates an intensive massage ideal for anti-cellulite treatments, e.g. on the thighs, buttocks and abdomen.
6. Cloud head - for massage of the back and other hard parts of the body.

WARNINGS

- When using the device for the first time or after a long break, be sure to charge the battery according to the instructions.
- If any abnormalities in the functioning of the body occur while using the device, stop using the device immediately and consult a doctor.
- If the product has been mechanically damaged, do not use it.
- Before charging, make sure that the charger and the charging socket of the product are dry.
- The massager cannot be soaked or immersed in liquid.
- The device is not a toy and should be kept out of reach of children.
- Packaging materials, such as plastic foil or other small parts, may be dangerous to small children.
- Do not use other massage devices at the same time while using this device.
- This product is not a medical device and cannot be used for therapeutic purposes.
- We recommend using the device for no longer than 15 minutes each time.
- The device contains a lithium-ion battery. It cannot be disassembled, modified, immersed in liquid or thrown into a fire, which may damage it or even cause the cell to explode.
- Do not use a power supply with a voltage and current other than those recommended by the manufacturer (5V, 1-2A) as this may damage the device. It is normal for the surface of the product to heat up during operation.



Important notes on safety of use

- Pregnant or recently pregnant women should use the massager after consulting a doctor
- People with hypersensitivity and skin diseases should consult a doctor about the possibility of using the device
- People with malignant tumors, heart disease, severe hypertension should be monitored by a doctor
- Do not use on children under 10 years of age.
- Disabled people require supervision when using this device
- Do not use in places with high humidity and high temperature, such as bathrooms, steam baths, saunas, drying rooms, etc.
- Do not use the product in the shower
- Do not use the product at temperatures below 0C
- For long hair, use a mesh cover to prevent it from getting tangled

Important. The device is equipped with a lithium-ion battery. It should be charged at least once every 3 months to maintain its full capacity. We recommend fully charging the battery before first use

Important. Due to the electronic components, the device must not be allowed to get wet (contact with water).

Battery charging

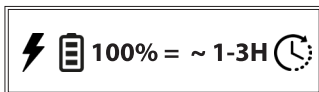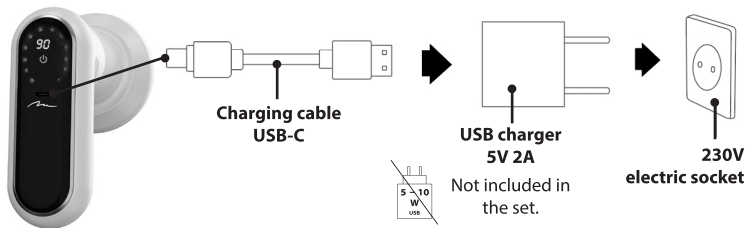Turning on, turning off and changing operating modes

- A short, single press of the touch button (3) displays information about the battery charge level in the range of 0-99.
- Pressing the touch button (3) for ~3s – turns the device on or off
- A short press of the button (3) changes the rotational speed of the massager (one of 6): L1, L2, L3, L4, L5, L6
- Successive presses of the button activate one of 3 vibration modes (at the highest speed): P1, P2, P3

Using the device and massage technique

The Anticell massager was created with the idea of effective reduction of any type of cellulite - regardless of its degree of advancement, location or origin (fatty, water, circulating cellulite). Daily, 10-15-minute treatments performed with the use of this massager are very effective: they break down accumulated fat cells, restore proper blood and lymph circulation.

The device emits infrared light during operation, gently heating the massaged area. It has a firming effect, improves skin elasticity and flexibility. It is effective in breaking down fat tissue and eliminating its excess. Reduces muscle tension, helps reduce spasms and swelling. It has an analgesic effect, relaxes and increases the feeling of relaxation of the body.

For the most effective massage treatment, the device should be held vertically just above the massage head. Do not hold it by the display part to avoid accidental touching of the button (3) and unwanted mode change.

Avoid applying too much pressure, which may cause pain or discomfort during the massage.

Optimal massage technique


Arms

Massage both arms evenly, spending the same amount of time on each. Massage the arms starting from the elbows towards the armpits, in small circles. This form of massage relieves muscle tension, shapes the body, firms and smoothes the skin.


Waist and belly

Massage the waist up and down. Then the waist from the outside to the inside. Perform a circular massage around the navel (clockwise). Then the abdomen, in large circles (also clockwise). This technique helps to effectively shape the waist and abdomen

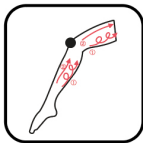Legs

Massage from the ankles upwards. Briefly massage the legs in straight lines upwards (always towards the heart). Then massage in small circles - from the bottom upwards (towards the heart). Special attention should be paid to those parts of the body where cellulite is present. Massaging in this way easily breaks down the fatty tissue, firms and smoothes the skin.

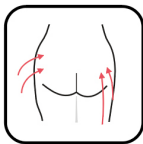Hip and buttock massage

Massage both parts of the body with movements from the bottom to the top. This type of massage is used to quickly reduce cellulite, lift the buttocks and firm the skin.

Other parts of the body

The massage should be performed in a circular motion or in a straight line - always towards the heart.

*Due to continuous development, product specifications and appearance are subject to change without notice.
For technical support, visit www.media-tech.eu.*

Specification:

- Wireless anti-cellulite massager
- Effective reduction of cellulite on the skin of the body and legs
- Improves the appearance and firmness of the skin
- Relaxation, stress relief, pain reduction, slimming
- 5 massage heads + fabric cover
- 6 speeds and 3 massage modes: 1800-3200 revolutions
- Gentle heating with infrared light
- LED display with touch control
- Power supply: USB-C 5V 2A socket (USB-C cable included)
- Battery: lithium-ion, up to 70 min of work on a single charge (in 10 min intervals)
- Charging time: ~ 3-4 hours
- Dimensions: 140x86x150mm

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**Environmental guidance:**

The device is marked with the symbol of a crossed out waste container, in accordance with the European Directive 2012/19/EC on waste electrical and electronic equipment. Products marked with this symbol should not be disposed of or disposed of after the end of their useful life together with other household waste. The user is obliged to dispose of waste electrical and electronic equipment by delivering it to a designated point where such hazardous waste is recycled. The collection of this type of waste in separate places and the proper process of their recovery contribute to the protection of natural resources. Proper recycling of waste electrical and electronic equipment has a beneficial effect on human health and the environment. In order to obtain information on where and how to dispose of waste electrical and electronic equipment in an environmentally safe way, the user should contact the relevant local authority, the waste collection point or the point of sale where he purchased the equipment.

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EN: The full text of the CE declaration of conformity is available at www.media-tech.eu.

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