# PL EMS NECK MASSAGER

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MT6534 V1.0

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
## User manual

MT6534 V1.0


*Congratulations on purchasing a high-quality EMS massager. Please read this user manual carefully to ensure full satisfaction and safety when using this device.*

**IMPORTANT!** The massager works based on EMS muscle electrostimulation technology with low-intensity current. People with implanted electronic circulatory support and heart function, as well as problems with arrhythmia or after a heart attack, should consult a doctor about the possibility of safe use of this device. People diagnosed with hypertension and atherosclerosis should consult a doctor before using the device.

## Appearance of the device
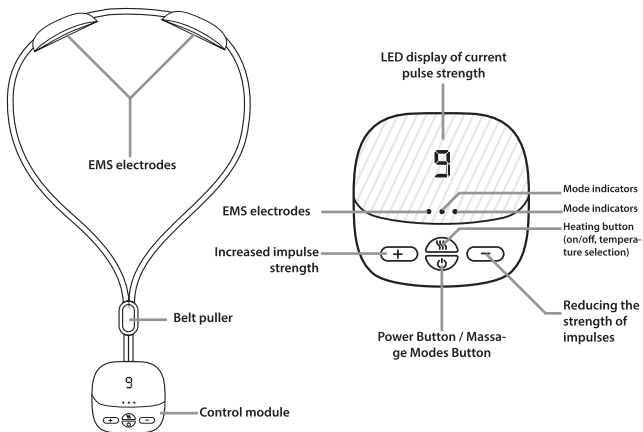
## Charging


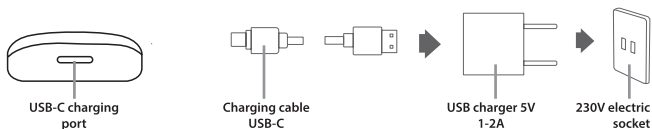
The device should be charged for about 2 hours. The LED screen will show the charge level until the letter F is displayed, which indicates that charging is complete.

**Important! Chargers with a voltage higher than 5V and a current higher than 2A cannot be used**

## How to put it on

### **IMPORTANT**

**Before using the massager, we recommend wetting the nape and neck with a damp towel and/or applying a small amount of water gel to the device's electrodes. This will prevent a tingling sensation of dry skin.**

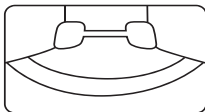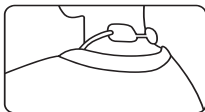Put on the massager so that the EMS electrodes adhere smoothly to the neck with their entire surface. For better effect and massage efficiency, use the belt tightener.

## How to use

### **IMPORTANT!**

**We strongly recommend starting with the lowest pulse intensity setting to accustom your muscles and skin to the effort. You should not exceed 2 sessions in one go (2 x 15 min).**

- Turn on the device with the power button by pressing and holding it for ~2s.
- Use the +/- button to select the pulse intensity (we recommend starting with level 1).
- Use the power button to select one of the 4 massage modes (short presses).
- The two mode indicator diodes will indicate the selected mode in the following order: green > green and red > red > 2 blue.
- The massage cycle lasts 15 minutes after which the device automatically turns off.
- To turn off, press and hold the power button for 2s.

## Warm compress

The device has a function of heating EMS electrodes. To turn it on, press the heating button. By pressing cyclically, set one of 3 temperatures (37-43°C) or turn off the heating.

You can turn on only the heating itself (warm compress), setting the EMS pulses to 0 with the +/- buttons

## WARNINGS

- When using the device for the first time or after a long break, be sure to charge the battery according to the instructions.
  - If any abnormalities in the functioning of the body occur while using the device, stop using the device immediately and consult a doctor.
  - If the product has been mechanically damaged, do not use it.
  - Before charging, make sure that the charger and the charging socket of the product are dry.
  - The massager cannot be soaked or immersed in liquid.
  - Do not use other massage devices at the same time while using this device.
  - This product is not a medical device and cannot be used for therapeutic purposes.
  - We recommend using the device for 15 minutes each time.
  - The device contains a lithium-ion battery. It cannot be disassembled, modified, immersed in liquid or thrown into a fire, which may damage it or even cause the cell to explode.
  - Do not use a power supply with a voltage and current other than those recommended by the manufacturer (5V, 1-2A) as this may damage the device.
  - It is normal for the surface of the product to heat up during use.
  - The device strap is covered with high-quality food-grade silicone that is neutral to the skin. However, people who are allergic to silicone should consult a doctor.
  - The EMS electrodes should be covered with water gel to avoid possible burns from too dry skin, especially at higher pulse strength settings.
- Driving while using an EMS massager is prohibited!**

## Important notes on safety of use

- Pregnant women should use the massager under the supervision of a doctor.
- Using the device after a stroke is prohibited.
- In the event of a fever above 38°C, the massager should not be used or the warm compress function should be turned off
- People with hypersensitivity and skin diseases should not use the device
- People with malignant tumors, heart disease, severe hypertension should be monitored by a doctor
- It is not recommended for use in children under 7 years of age.
- Disabled people require supervision when using this device
- Do not use in places with high humidity and high temperature, such as bathrooms, steam saunas, saunas, drying rooms, etc.
- Do not use the product in the shower
- Do not use the product at temperatures below 0°C
- Do not use the product near devices emitting microwaves (min. 1 meter)
- Do not short-circuit the EMS electrodes during device activity

**Important. The device is equipped with a lithium-ion battery. It should be charged at least once every 3 months to maintain its full efficiency. We recommend fully charging the battery before first use.**

**Important. Due to the electronic components, the device cannot be allowed to get wet (contact with water).**

*Due to continuous development, product specifications and appearance are subject to change without notice. For technical support, visit [www.media-tech.eu](http://www.media-tech.eu).*

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**EN Environmental guidance:**  
The device is marked with the symbol of a crossed out waste container, in accordance with the European Directive 2012/19/EC on waste electrical and electronic equipment. Products marked with this symbol should not be disposed of or disposed of after the end of their useful life together with other household waste. The user is obliged to dispose of waste electrical and electronic equipment by delivering it to a designated point where such hazardous waste is recycled. The collection of this type of waste in separate places and the proper process of their recovery contribute to the protection of natural resources. Proper recycling of waste electrical and electronic equipment has a beneficial effect on human health and the environment. In order to obtain information on where and how to dispose of waste electrical and electronic equipment in an environmentally safe way, the user should contact the relevant local authority, the waste collection point or the point of sale where he purchased the equipment.

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## FEATURES

- EMS electrostimulation neck muscle massage device
- Improve blood circulation, reduce pain, neck stiffness and tense muscles
- Warm compress function –37–43°C (3 levels)
- Very easy operation and design for use anytime and anywhere
- Dual pulse technology
- 4 massage modes
- 9 levels of intensity
- Levels 1-3: normal relaxation massage
- Levels 4-6: reduction of chronic cervical spine pain
- Levels 7-9: reducing the inconvenience of "uncomfortable pillow"
- Control module with LED display
- Pulse power: ~5W
- Power supply: rechargeable lithium-polymer battery
- Operation: up to 150min (in 15min intervals) on a full charge
- Battery charging: 5V 1A USB-C cable ~2 hours
- Dimensions: 4,7x4,7x1,6cm (controller), max neck circumference ~78cm
- Weight: 0.57kg