# EMS FOOT MASSAGER

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
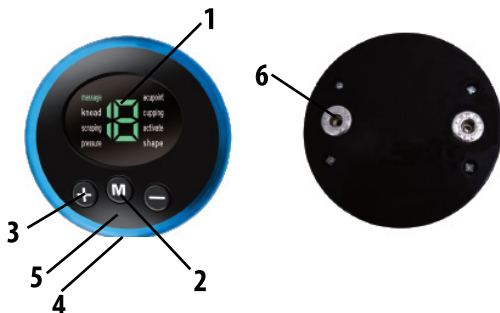
## Mat and remote control

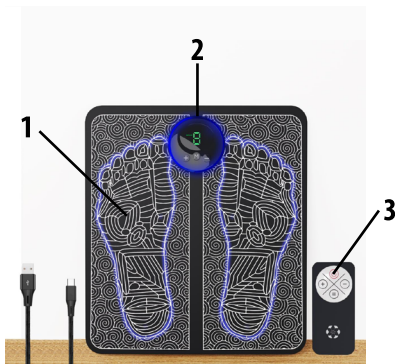1. The location of the user's feet so that they adhere to the largest possible surface
2. Mat power latches for controller placement
3. Remote control with +/- adjustment buttons, ON/OFF button and mode change button

## Charging and installation

Before starting the massage, we recommend charging the controller battery to 100%. To do this, connect the socket [4] with the included USB/microUSB cable to the 5V/1-2A charger or USB port on the computer. Do not use the so-called fast chargers. **Do not use the device while charging**  
**The charging of the device is signaled by the LED [5], which will turn off after charging is complete.**

Remember to **fully** charge the device after each use and charge it at least once every 3 months, even when the device is not in use, to keep the battery in good condition.

Remove the battery discharge protection from the remote control. Connect the controller to the charging latches so that a "click" is heard

## Usage

- a. Moisten the soles of your feet with a damp towel for a better massage effect. The procedure should be performed sitting down so that the leg muscles are unburdened and relaxed. Do not put your full weight on the mat
- b. Place your feet on the designated places on the mat so that they adhere as much as possible
- c. Use the [+] button to turn on the device
- d. Use the mode button [M] to select one of the massage modes
- e. Use the +/- buttons to set one of the 19 levels of intensity.

**The mat will turn off automatically when one foot is removed.**

The massage involves the use of mild electrical impulses (EMS) to stimulate the muscles of the feet and calves. This helps to improve blood circulation, relieve fatigue and muscle tension, and stimulate the acupressure points on the feet. The mat has different massage zones, covering the sole, heels, toes and arches of the feet.

The first use may require getting used to a slight sensory discomfort, especially in people with sensitivity to the so-called. "ticklers".

We recommend that you select the massage intensity level "1" the first time and find the treatment modes that best suit you.

EMS treatment is recommended maximum once a day, up to 15 minutes of continuous massage. The device has a protection that automatically turns off the treatment after about 15 minutes.

**IMPORTANT!**

**The MT6527 is not a medical device and cannot replace recommended therapy in medical centers. In case of doubt, it is recommended to consult a doctor.**

**IMPORTANT!**

**The mat works based on EMS muscle electrostimulation technology with low current. People who have implanted electronic circulatory support and heart function, as well as problems with arrhythmia or after a heart attack, should consult a doctor about the safe use of this device.**

**Before using the device, medical consultations should be performed by people with diagnosed arterial hypertension and atherosclerosis.**

*Due to continuous development, product specifications and appearance are subject to change without prior notice. For technical support, please visit [www.media-tech.eu](http://www.media-tech.eu).*

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**EN Environmental guidance:**

The device is marked with the symbol of a crossed out waste container, in accordance with the European Directive 2012/19/EC on waste electrical and electronic equipment. Products marked with this symbol should not be disposed of or disposed of after the end of their useful life together with other household waste. The user is obliged to dispose of waste electrical and electronic equipment by delivering it to a designated point where such hazardous waste is recycled. The collection of this type of waste in separate places and the proper process of their recovery contribute to the protection of natural resources. Proper recycling of waste electrical and electronic equipment has a beneficial effect on human health and the environment. In order to obtain information on where and how to dispose of waste electrical and electronic equipment in an environmentally safe way, the user should contact the relevant local authority, the waste collection point or the point of sale where he purchased the equipment.

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