


Quick Guide

***Statement**
Compared with manual toothbrush, this sonic electric toothbrush has as much as 7 times of plaque removal effect. When the first time use, tickling or slight bleeding may occur. It is because the gums are not used to this kind of efficient clean. But this phenomenon usually will disappear within 1 to 2 weeks. If this phenomenon lasts for more than 2 weeks, it indicates your gums are abnormal. In this case, please stop using this sonic electric toothbrush and consult the dentist.

*Warning

1. Please use under guidance of the dentist if users had dental correction, dental repair, dental implant.
2. This product is designed to clean teeth, gums and tongues. Please don't use for any other purposes. If feel any discomfort or pain, please stop using and see a doctor.
3. Don't use a deformed or bristle- bent toothbrush head. Advise to change the toothbrush head every 3 months. If there are signs of wear, the toothbrush head should be replaced in advance.
4. Do not knock or drop the product. If it is damaged due to a collision or falling, please stop using it.
5. Do not use this sonic electric toothbrush in a bath or shower.
6. This product contains a built- in rechargeable battery. Please do not put it into the fire, and don't charge, use or place it at high temperature.

Product Specifications

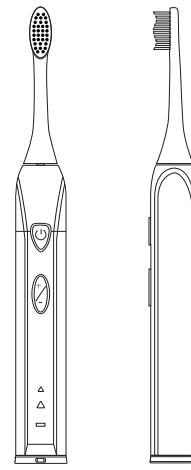


Product Size
246mm*33mm



Input Voltage
DC5V


Accessories

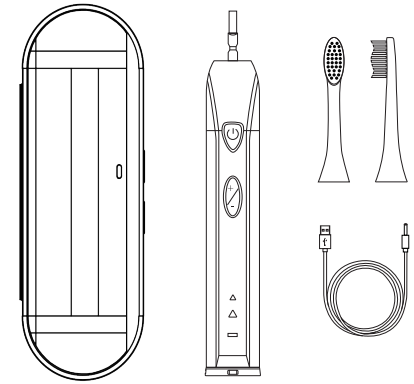

| Storage box | Device body | Brush head | USB cable |
|-------------|-------------|------------|-----------|
| 1 | 1 | 2 | 1 |

Charging

- Open the bottom waterproof cover, and use USB cable to charge.

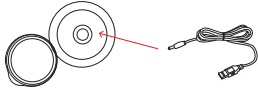

- The power light is breathing red when charging. It will turn green after fully charging.
- When low power, The power light will be red and quick flashing. Fully charging need 6 hours(the first fully charging need 12 hours).
- Suitable for DC5V input (round jack).
- Before charging, please make sure that the charging socket and the power plug are completely dry, then connect to the power supply.

Precautions

Prohibited items

- Please keep away from strong electromagnetic fields when charging, such as induction cookers.
- Do not charge in the bathroom. Do not place the charger above or near the water containers, such as tub, wash-basin, pool and so on.
- Do not place metal objects, such as coins or paper clips, on the charger.
- If will not use it for a long time, please cut off the power supply to avoid causing fire or danger.

Function Description

Function of the keys


- On/Off Key: (2 functions)

1 Starting up: short press.

2 Suspend: short press On/Off Key. Automatically shutdown after 5 seconds.

- Speed key:

Total 2 speeds, +/- control the speeds. (It can adjust the speeds under on or off status.)

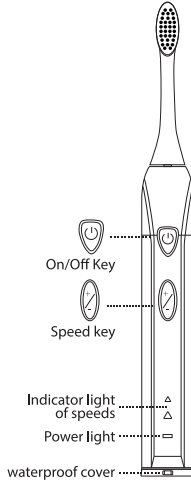

2 minutes smart timing

- Every 30 seconds to pause for 1 second to remind reversing brushing the teeth. Brushing time is 2 minutes.

Maintenance

Clean brush head:

- Please clean brush head after each use, and let it naturally dry.
- Do not pull the bristles, which may cause the bristles deformation and falling off.


Clean the device body:

- Scrub with a soft cloth and let it dry naturally.
- Do not immerse it in water.
- Don't use higher than 40℃ water or acid cleaning fluid to clean.

Replace the brush head:


- The dentist recommends replacing the brush head every three months.
- Please use the manufacturer's recommended toothbrush head.

Brushing time

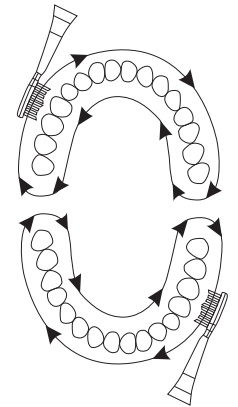


Dentists recommend that brush your teeth for 2 minutes each time. During the brushing, every 30 seconds to have one second pause to remind users to change teeth area to brush.


The 1st step: Place the brush head toward the front, and then press down the brush head to the end until only a gap about 1 mm between the brush and the handle.

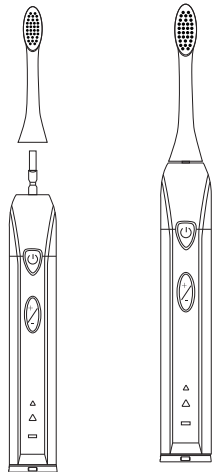


The 2nd step: Squeeze the appropriate amount of toothpaste on the toothbrush, place the toothbrush gently at the 45° junction of the teeth and gums.

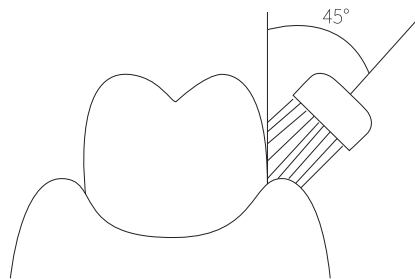


The 3rd step: Press the On/Off key to start cleaning your teeth.

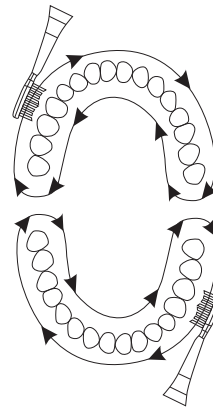

The 4th step: To prevent toothpaste foam splash, please stop the toothbrush brushing before taking it out from mouth.

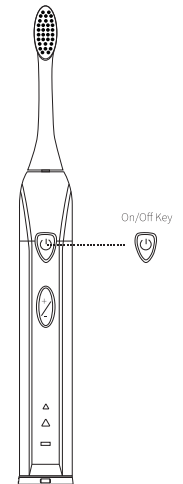

The 5th step: After finish brushing, rinse the toothbrush head with flowing clean water. Shake water off the brush head, and wipe clean the toothbrush handle with a soft cloth.

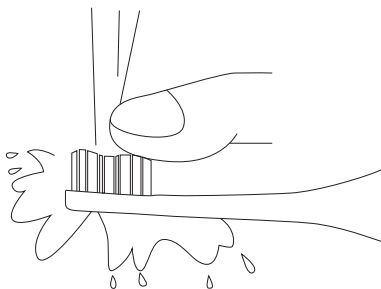

The 6th step: Place the sonic toothbrush vertically.

